


PRISMATIC DENTIN



PRISMATIC FLUORESCENT

ENAMEL MASTER

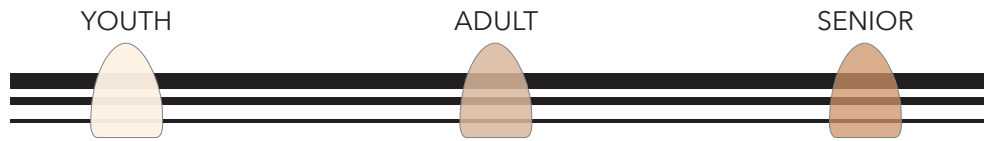

RESTORATION AGE	PART OF	PRISMATIC DENTIN	PART OF	SISTEM INTERNAL DENTIN	CHROMA	PART OF	OTHER MASSES	RESULT
TEEN	1	YOUTH	1	SID LIGHT	A	2	Transpa NEUTRAL	A1
					B			B1
					C			C1
					D			D2
YOUNG	1	YOUTH	2	SID LIGHT	A		A2	
					B	B2		
					C	C2		
MATURE	1	ADULT	2	SID DARK	A		A3,5	
					A	A3,5		
					B	B3		
					C	C3		
					D	D3		
SENIOR	1	SENIOR	2	SID DARK	A		A4	
					B	B4		
					C	C4		
					D	D4		



By *Nunzio Rosa*

1	1	PARTS OF	FLUO Safari 607	+	2	PARTS OF	FLUO Sand 622	=	brighten and vivid colors
2	2		FLUO Orange 623	+	1		FLUO Red Passion 616	=	abrasion, cuspid internal
3	1		FLUO Chocolate 106	+	1		FLAT Brown 508	=	incisal areas, rift zones, heavy colored interproximal areas
4	1		FLUO Atlantic Blue 621	+	2		FLUO London 617	=	deep blue, great as incisal effect
5	1		FLUO Atlantic Blue 621	+	2		FLUO Blue 619	=	lower the value
6	2		FLUO Arancio 623	+	1		FLUO Chocolate 109	=	rust color, perfect for darker necks and roots
7	3		FLUO Giallo 601	+	1		FLUO Chocolate 109	=	ochre
8	1		FLUO Atlantic Blue 621	+	1		FLUO Blue 619	+	1 PARTS OF FLUO Meteorite 612 = deep dark gray, tetracycline effect

